

A clear, multi-lumen medical infusion set with blue tubing, red and green connectors, and blue clamps, arranged in a loop.

Mediplus TIVA Sets

Dedicated infusion delivery systems

40

Years of innovative patient care

Mediplus TIVA Sets

Dedicated infusion delivery systems

Mediplus TIVA Sets were designed in conjunction with anaesthetists to meet the highest TIVA standards as outlined in the SIVA, Euro SIVA and AAGBI guidelines. These high-performance sets offer a wide range of benefits including:

Enhanced performance

- Using a dedicated line reduces the risk of leakage associated with multiple connections
- Integrated anti-siphon valve and high flow anti-reflux valve prevent free flow and back flow
- Wide bore IV gravity line enables high fluid flow to ensure continuous infusion
- Minimal dead space designed to help prevent potential overdose in cases of accidental bolus
- Micro bore tubing to prevent kinking and occlusion
- Siamese tubing on double and triple infusion lines to prevent entanglement

Greater patient comfort

- The soft distal end minimises stress at the cannula site to help avoid displacement

Cost savings

- Micro bore tubing minimises post-operative drug wastage

Wide range

- Single, double, triple, quadruple and quintuple infusion lines available in a range of lengths

“We used the Mediplus TIVA giving set for over 15 years without incident or failure. As a cost saving measure we changed to an alternative giving set, but within four months we found that product to be inferior and we have now returned to our tried and trusted product”

Bruce Rodda, Consultant Anaesthetist
Mid Essex Hospital Services NHS Trust



Mediplus TIVA Sets

Dedicated infusion delivery systems


Key features



Anti-Siphon valves to prevent back flow and free-flow



High flow anti-reflux valve to prevent back flow

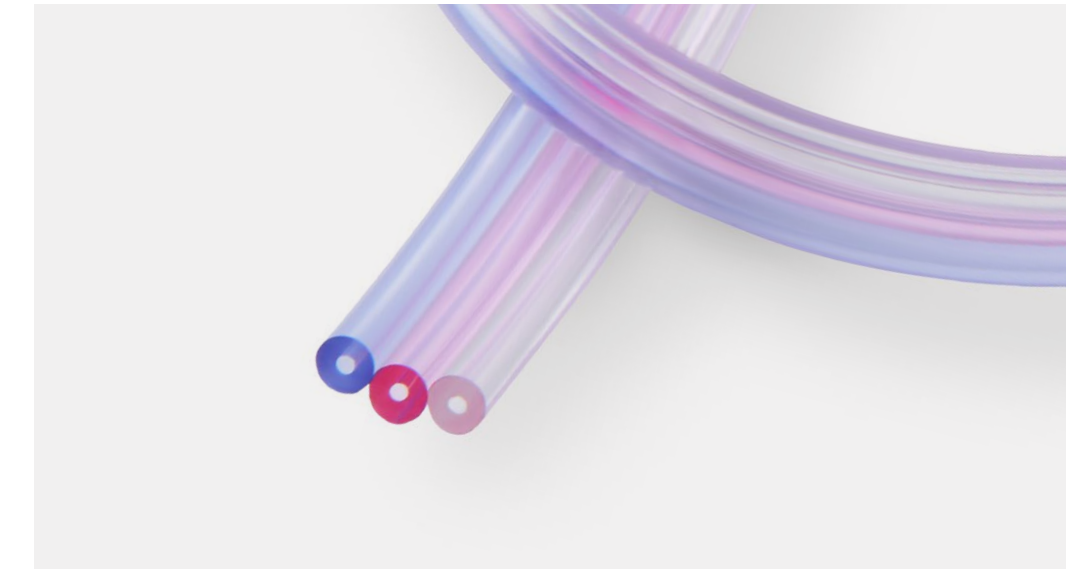


Micro bore siamese tubing for a low priming volume but large external diameter to prevent kinking/occlusion



Red end caps always highly visible



Soft distal end to provide minimal stress at the cannula site

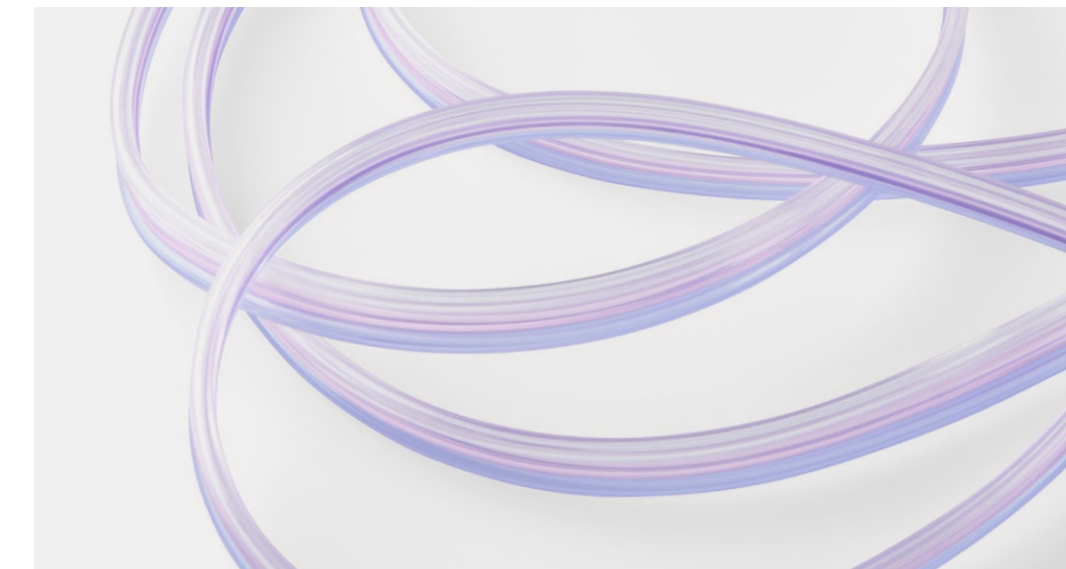

Siamese configuration to prevent entanglement of tubing



Dual TIVA Set



2-way TIVA Set



3-way TIVA Set



4-way TIVA Set



5-way TIVA Set



Triple Lumen Extension Set

Ordering information

Product Code	Product Description	Length	Box Quantity
6611	Single lumen extension set: 1 x anti-siphon valve (200cm/78")	2.0m	50
6616	Single lumen extension set: 1 x anti-siphon valve (600cm/236")	6.0m	25
6610	Double lumen extension set: 2 x anti-siphon valve (200cm/78")	2.0m	50
6604	Double lumen extension set: 1 x anti-siphon, 1 x high flow back check valve (200cm /78")	2.0m	50
6603	Triple lumen extension set: 2 x anti-siphon, 1 x high flow back check valve (200cm/78")	2.0m	50
6605	Triple lumen extension set: 2 x anti-siphon, 1 x high flow back check valve (250cm/98")	2.5m	25
6606	Triple lumen extension set: 2 x anti-siphon, 1 x high flow back check valve (300cm/118")	3.0m	25
6607	Triple lumen extension set: 2 x anti-siphon, 1 x high flow back check valve (600cm/236")	6.0m	10
6633	Triple lumen extension set: 3 x anti-siphon valve (600cm/236")	6.0m	10
6630	Quadruple lumen extension set: 3 x anti-siphon, 1 x high flow back check valve (220cm/87")	2.2m	25
6631	Quadruple lumen extension set: 3 x anti-siphon, 1 x high flow back check valve (300cm/118")	3.0m	25
6635	Quintuple lumen extension set: 4 x anti-siphon, 1 x high flow back check valve (250cm/98")	2.5m	25
6636	Quintuple lumen extension set: 4 x anti-siphon, 1 x high flow back check valve (300cm/118")	3.0m	25

Clinical data

Guidelines for the safe practice of total intravenous anaesthesia (TIVA) Joint Guidelines from the Association of Anaesthetists and the Society for Intravenous Anaesthesia A. F. Nimmo, A. R. Absalom, O. Bagshaw, A. Biswas, T. M. Cook, A. Costello, S. Grimes, D. Mulvey, S. Shinde, T. Whitehouse and M. D. Wiles.
Tips and Tricks to optimize Total Intravenous Anesthesia L. Barvais, F. A. Lobo, F. H. M. Engbers, M. G. Irwin, T. W. Schnider and S. Schraag.



Mediplus are experts in the delivery of high performance medical devices that support positive patient outcomes in the following fields -



anaesthesia



gynaecology



operating room



urology

Mediplus Ltd.
+44 (0) 1494 551200
help@mediplusuk.com
The Gateway Centre, Coronation Road,
High Wycombe HP12 3SU
www.mediplusuk.com

INVESTORS IN PEOPLE®
We invest in people Gold

